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ENGLISH

You have purchased a RHINO product, which guarantees quality and durability, since the BABAIN-180NG scale was manufactured according to American quality standards. Keep this instruction manual, it will help you understand the operation of the RHINO BABAIN-180NG scale will be useful for future reference.

TECHNICAL SPECIFICATIONS

ACCURACY

Medium III

CAPACITY / DIVISION

400 lb / 0.1 lb

POWER

3 Vcc
(2 batteries AAA 1,5 Vcc)

TEMPERATURE

32 °F - 104 °F

QUICK GUIDE TO USE

- Remove the scale from its packaging, make sure it's in perfect conditions and that includes all its accessories (manual and 2 AAA batteries).
- Place the scale on stable and leveled surface, be sure to place batteries correctly any movement no matter how small can affect accuracy.
- Place your feet in the center of the scale.

CAUTIONS

In case of need of technical assistance, please contact one of the authorized Rhino service centers.

- The display will show "**Lo**" when it is necessary to replace the batteries.
- The display will show "**Err**" when the scale is overloaded. Remove excess weight to avoid damage.
- Place the scale in a cool, dry place.
- For domestic use only.
- DO NOT USE for commercial purposes.

WARNINGS

In case of need of technical assistance, please contact one of the authorized Rhino service centers.

Do not drop the scale suddenly.

To clean the scale use a damp cloth or sponge, never a direct stream of water.

Protect the scale from insects, do not spray, only solid products against insects.

Do not open the scale yourself or the warranty will be void.

Do not direct air flow from fans or heaters directly to the scale during use.


Do not expose your scale to dripping or splashing liquids.

BATTERY INSTALATION

- Open the lid underneath the scale as shown in the diagram.
- Install the batteries taking into account the polarity indicated in the compartment.
- Replace the battery compartment cover.

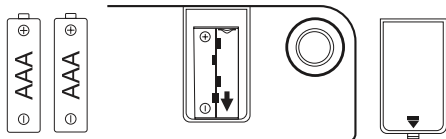WEIGHING PROCESS

- To ensure the accuracy of the scale, place it on a hard, flat surface.

NOTE. Avoid placing it on a slippery surface.

- Gently step on the scale, it will turn on automatically.
- Stay on the scale without move until the weight appears in the display, when the numbers on the display flash 3 times and then remain fixed, it will indicate its weight.

WARNING TO USE

- DO NOT use with medical implants as pacemakers.
- DO NOT stand on the edge of the scale or jump on it
- DO NOT overload the scale (MAX. 180kg / 396lb).
- DO NOT drop the scale or objects on it, as this could damage the sensors.
- NOT recommended for women pregnant.
- Any measurement obtained with this device is for reference only and should not be considered as a medical opinion.
- Always place the scale on a hard, dry and flat surface before of measurement
- Make sure your feet are dry before stepping on the scale.
- Remove batteries if scale is not used for an extended period of time.
- The user under 6 years old, just use weight function without BMI (body mass index) and body fat, data.
- The user under 10 years old, use only weight and BMI function, no body fat, data etc.

INDICATION SYMBOLS

Lo - Low battery, change batteries.

Err - Overload.

C - An error has occurred during measurement. Get off the scale and repeat the process.

DOWNLOAD THE APP

To use **FITDAYS** app you must to download the app to your smartphone. You can do it from Google Play or App Store.


NOTE. The app is compatible only with IOS 8.0 or higher Android 5.0 or higher operating systems and only supports one smartphone device at a time.

- After downloading the app, check that your bluetooth is turned on.
- Enter the app from your phone or tablet and register your email.
- Read and accept the privacy policies.
- Add personal data like sex, age, height etc (Enter incorrect data may result in inaccurate measurements).

UNIT CONVERSION

The unit change button **UNIT** is located underneath the scale, as shown in the diagram. With the button **UNIT** you can change the unit of measurement, kilograms (kg) to pounds (lb) and vice versa.



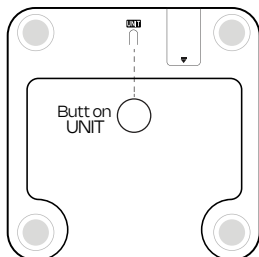
HOW TO USE THE APP

- Kick the scale to activate the display. Wait until the display shows **"0.0"** (kg / lb).
- With your app open on the home screen, place your bare feet on top of the metal bars.

NOTE. You must place your bare feet on the metal bars.

If you step on the scale with socks or shoes, only the weight and BMI (body mass index) will be displayed on the main page.

- The numbers in the screen will be blinking while the scale is analyzing date.
- Keep your feet on the scale till the weight showed in the app be stable then the scale will show the obtained information.
- You will find the result such as BMI, visceral fat, body water, skeletal muscle, muscle mass, body age, protein, subcutaneous fat, muscle frequency, fat loss and BMR. Also you can see the weekly, monthly and year result.

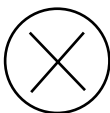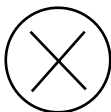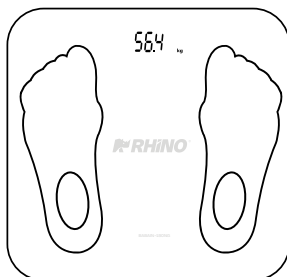USER



- Add or change users by clicking on **"MY PROFILE"** or click on **"NEW USER"**.
- The scale and the app recognize the variation of results, and gives you the option to change the user if the data is different from the selected user.

NOTE. Each account can support 24 users.

BABY WEIGHING MODE

- Click **"OPTIONS"** then select **"BABY MODE"** in the app.
- Get on the scale so the app is determined its weight before carrying the baby.
- After register your weight, hold the baby and wait stabilize the scale to get the baby's weight on the application.

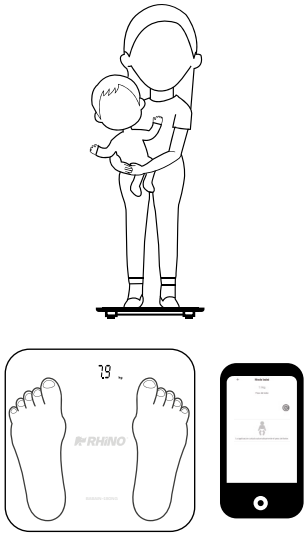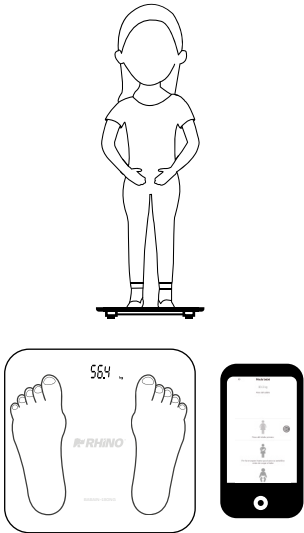TEAM BODY PROBLEMS AND SOLUTIONS

PROBLEM	SOLUTIONS
Bluetooth not connected	The scale already has bluetooth activated once turned on. Check the IOS version of your smartphone must be IOS higher than 8.0 or higher / Android 5.0 or higher. Restart your mobile device. Remove the batteries from the scale and put them back. NOTE. The scale only supports one smartphone device at a time.
Body fat data is not measured	User must be barefoot to measure body fat. Make sure the scale and your feet are dry.
The scale is producing inconsistent weights	Check if the scale is on a hard flat surface. Check each foot of the sensor to make sure nothing is stuck to the bottom of the sensor.
The app shows the wrong ratio of muscle	The body parameter setting is incorrect, check to make sure I correctly enter the sex, age and height data.
The scale does not turn on	Remove the batteries and put them back. Be sure to insert the batteries into the correct address. Carry out the battery change.

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WARRANTY

Rhino Maquinaria warrants to customers who purchase any Rhino products from an authorized seller that Rhino products will be free of defects in workmanship for the duration of 2 years from the date of purchase.

During the 2 years warranty period, Rhino Maquinaria will replace the product that has a manufacturing defect as long as it has been determined by Rhino Maquinaria.

This warranty does not apply to defects caused by misuse or alterations of the product, nor does it apply if the damage was caused by repairs or service provided or attempted by anyone other than Rhino Maquinaria.

To make this warranty effective, you must first contact our Customer Service by dialing (619) 393-1079 or by mail to info@rhinoglobal.net, once you have contacted our Customer Service they will request from you a warranty claim by email.

Rhino Maquinaria will review the warranty claim. If the warranty claim is found to be valid, please return the product with a copy of the invoice and order ID to the of the invoice and order ID to the indicate address.

Upon receiving the product, Rhino Maquinaria will sent a replacement product to customer's address.

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Rhino Maquinaria S.A. de C.V.